# HC-DT11-12

**MANUAL (p. 2)**  
Digital Thermometer

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
## ENGLISH

### What is a “normal” temperature?

Although the generally accepted “normal” temperature reading is 37.0°C/98.6°F, temperature readings can vary from 36.1°C/97.0°F to 37.2°C/99.0°F and still may be considered “normal”. Variations in temperature can be attributed to activities such as exercise, smoking, eating and drinking. Your temperature is lower in the morning than in the afternoon. Other variations may be due to the location of the temperature reading. While the oral temperature follows the guidelines previously described, a rectal temperature is generally 0.5°C /1.0°F higher. Conversely, an axillary (under the arm) temperature will be 0.5°C /1.0°F lower.


### Specifications

|                                          |                                                                                   |
|------------------------------------------|-----------------------------------------------------------------------------------|
| Measurement range                        | 32.0°C (89.6°F) ~ 42.9°C (109.2°F)                                                |
| Low temperature display                  | Temperature <32.0°C (89.6°F): display: Lo°C (Lo°F)                                |
| High temperature display                 | Temperature >42.9°C (109.2°F): display: Hi°C (Hi°F)                               |
| Display resolution                       | 0.1°C (32.2°F)                                                                    |
| Accuracy                                 | 35.0°C (95.0°F) ~ 39.0°C (102.2°F): ±0.1°C (0.2°F)<br>Other range: ±0.2°C (0.3°F) |
| Fever beep Sound<br>[>37.8°C (100.0°F)]  | 30 short beep sounds in 10 seconds                                                |
| Normal beep sound<br>[≤37.8°C (100.0°F)] | 10 long beep sounds in 10 seconds                                                 |
| Auto power off                           | 8 ± 2 minutes                                                                     |
| Battery                                  | LR41 or SR41                                                                      |
| Life of battery                          | 2 years                                                                           |
| Operating conditions                     | Temperature: 10°C (50°F) ~ 40°C (104°F)<br>Humidity: 10%~95% RH                   |
| Storage conditions                       | Temperature: -20°C (-68°F) ~ 60°C (140°F)<br>Humidity: 10%~95% RH                 |
| Weight                                   | HC-DT11-12: 11g    Frog: 16 g    Duck: 14 g                                       |
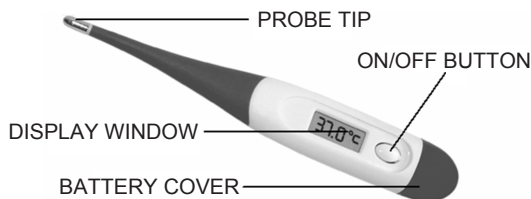
### How to take your temperature

1. Press ON/OFF button to activate. The unit will beep and display Lo. This LCD display test will last for about 2 seconds.
2. Then, the temperature taken last time will be displayed with an “M” mark in the right corner.
3. 37.0°C (or 98.6°F) will appear.
4. When the thermometer displays Lo and flashing °C (°F) it's ready for temperature measuring. If the room temperature is higher than 32.0°C (89.6°F), the room temperature will be displayed instead of Lo.
5. The beep will sound when temperature measuring is complete.
6. Press the ON/OFF button to turn off the unit.  
(The unit will automatically power off approximately ten minutes after use.)

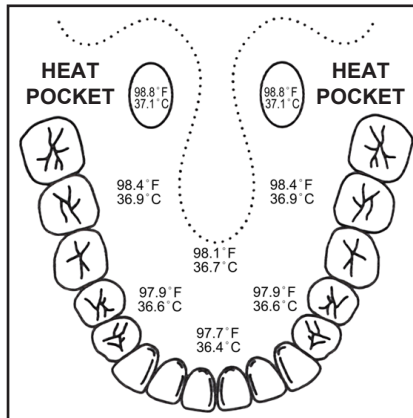Located on the bottom (floor) of the mouth at the base of the tongue

### Taking an oral temperature

Place the probe tip under the tongue as near as possible to the heat pocket.

Approximate measuring time: 30 seconds.

Note: Please remember that the oral temperature is a little higher than external temperature (about 0.3°C or 32.54°F).

Suggestion: Close your mouth at least 2 minutes before taking the temperature for a more accurate measuring.

### Taking an underarm temperature

Wipe the underarm with a dry towel and place the probe tip under the arm so the tip is touching the skin with the thermometer perpendicular to the body. Put your arm across the chest so the probe tip is covered properly.

Approximate measuring time: 10 seconds.

Note: Underarm temperatures are generally 0.5°C or 32.9°F lower than oral readings.

Suggestion: Close your armpit at least 5 minutes before temperature taking for a more accurate measurement.

### Taking a rectal temperature

Apply a water-soluble lubricant to the probe cover. Gently insert the probe (MAXIMUM ½") into the rectum.


Approximate measuring time: 20 seconds.

Note: Please remember that the rectal temperature is a little higher than external temperature (about 0.5°C or 32.9°F).

### Replacing the battery

The battery in your digital thermometer needs to be replaced when "⎓" appears in the lower right corner of the display. Replace the battery with a 1.55 V, LR41 alkaline or SR41 silver oxide type or equivalent.

### To replace the battery

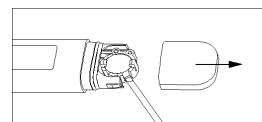1. Pull the battery cover away from the unit.
2. Remove the battery holder.
3. Use a non-metal pointed tool to remove the battery and insert the new battery with the positive (+) side facing up.
4. Replace the battery holder.
5. Carefully slide the battery cover back into place.

Warning: Please properly replace the battery away from small children and heat.

### Cleaning instructions

Before and after each use, clean the thermometer with a soft cloth and isopropyl alcohol diluted with water, or wash with COLD soapy water. The unit can be submerged in water or alcohol (waterproof model only). Do not soak for an extended period time or sterilize by boiling method, gas or steam autoclaves.

| Solution          | Dilution with water | Method                               |
|-------------------|---------------------|--------------------------------------|
| Ethanol           | 70-80%              | Wipe the probe with the soaked cloth |
| Isopropyl Alcohol | 70-90%              | As above                             |

### Precautions

1. The thermometer should only be used under the supervision of an adult. It is not a toy for children. Please keep it out of reach of children.
2. Do not walk, run or talk while measuring the temperature.
3. Store the unit in the protective case when not in use.
4. Do not store the unit where it is exposed to direct sunlight, dust or humidity. Avoid extreme temperatures.
5. Dropping or subjecting your thermometer to strong shocks should be avoided.
6. Do not attempt to disassemble the unit, except to replace the battery.

### Safety precautions:

**CAUTION**  
RISK OF ELECTRIC SHOCK  
DO NOT OPEN

To reduce risk of electric shock, this product should ONLY be opened by an authorized technician when service is required. Disconnect the product from mains and other equipment if a problem should occur. Do not expose the product to water or moisture.

### Maintenance:

Clean only with a dry cloth. Do not use cleaning solvents or abrasives.

### Warranty:

No guarantee or liability can be accepted for any changes and modifications of the product or damage caused due to incorrect use of this product.

### General:

Designs and specifications are subject to change without notice.

All logos brands and product names are trademarks or registered trademarks of their respective holders and are hereby recognized as such.

Keep this manual and packaging for future reference.

### Attention:

This product is marked with this symbol. It means that used electrical and electronic products should not be mixed with general household waste. There is a separate collections system for these products.

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Copyright ©
**NEDIS B.V.**